RELAXXNOW

Shiatsu + Tapping

RXM300

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Neck and Shoulder Massage Belt

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DEUTSCH

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ENGLISH

Technical Data	14
Safety and maintenance issues	14
Product features	17
Method and function of buttons	18
Cleaning & maintenance	19
Troubleshooting	20
Disposal	21

Thanks for purchasing this product. Please read this manual instruction carefully in order to know how to use this product correctly before operation. Please pay more attention to some important safety warning and keep this manual book properly for future use.

Note: Our Company reserves the right to revise the design and description of this product without any further notice. Please refer to the actual color of this product.

TECHNICAL DATA

Rated Power: 85 W

Net Weight: 2.8 kg

Input Voltage 220V-240V/ AC 60Hz

Rated working time: 15 minutes

Scope of delivery: Massage device | AC adapter | Manual

How to use: Insert the power plug into the socket

Product dimensions: 128 cm x17 cm x 16,5 cm

SAFETY AND MAINTENANCE ISSUES

1. Important safety warning:

- Do not allow any foreign objects to enter the space between the massage head and the casing.
- Do not allow the product to come into contact with water to avoid any possible accidents that may result from a short circuit.
- If the main cable or plug is damaged, do not disassemble and attempt to repair it by yourself. Please contact your distributor.
- Ensure that the main power is suitable for this product.
- Please disconnect the product from the main supply after use or before cleaning to avoid damage to the equipment or personal injury.
- Please use this product according to the user manual.
- Do not use any spare parts and attachments other than those recommended.
- Do not use this product if its casing or main cable is damaged.
- Do not use this product if you are under the influence of alcohol or feeling unwell.
- Do not use this product with too much force to avoid injury.
- Do not keep this product on for more than 15 minutes to avoid overheating.
- Do not use this product within one hour of eating.
- Please use this product carefully according to your doctor's advice if you have a neck or back injury or if you have had an operation recently.

2. Working environment:


- Do not use this product in a highly humid environment such as a bathroom.
- Do not use this product immediately after there has been a significant change in the environmental temperature.
- Do not use this product in an environment where there are lots of dust or corrosive gases.
- Do not use this product when it is covered by a blanket or pillow, etc.

3. People that should not use this product:

- People with contagious skin conditions.
- People with osteoporosis.
- People with heart diseases and with electronic medical instruments like a heart pacemaker.
- People with high fever.
- Women that are pregnant or menstruating.
- Those who are injured or suffering from any body surface disease.
- The appliance has a heated surface. People insensitive to heat must be careful when using the appliance.
- Children under 14-years or people mentally unstable without being supervised are not allowed to use this machine.
- Those who need rest according to a doctor's advice or those who do not feel well.
- Do not use this product when you are wet.
- Do not use this product on a pet.

4. Safety precautions:

- Do not plug or unplug this product when your hands are wet.
- Do not pull on the cable when unplugging. Do not plug or unplug carelessly.
- This appliance is not intended for use by people (including children) with reduced physical, sensory or mental capabilities or lack of experience and knowledge, unless supervised during use by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- Do not damage the cable or change the circuit of the product.
- Please avoid using this product in the case of a power outage in order to avoid any body injury upon sudden power resumption.
- Please stop using this product and consult your local supplier if you notice anything abnormal with the product during use.
- Please stop using this product and consult your doctor for advice.
- If the supply cord is damaged, it must be replaced by the manufacturer, service agent or similarly qualified person in order to avoid a hazard.
- The appliance has a heated surface. People insensitive to heat must be careful when using the appliance.
- Please do not pull controller connecting line, in case the power cord is fractured or when plug has bad contact. When pulling out the power plug, please hold the plug pulled back, rather than directly drawing the power cord.
- When the whole machine stops working, please immediately stop using it and pull out the power plug.
- This product is prohibited in the following locations: bathroom and kitchen, high temperature lights, damp, high voltage and the source of fire and water.
- Keep far away from sharp objects, the machine is not be used when its encloser is damaged.
- Do not force it with fierce pressure or put body pressure on the massage mat. Don't place hard objects on it; don't load or tread it.

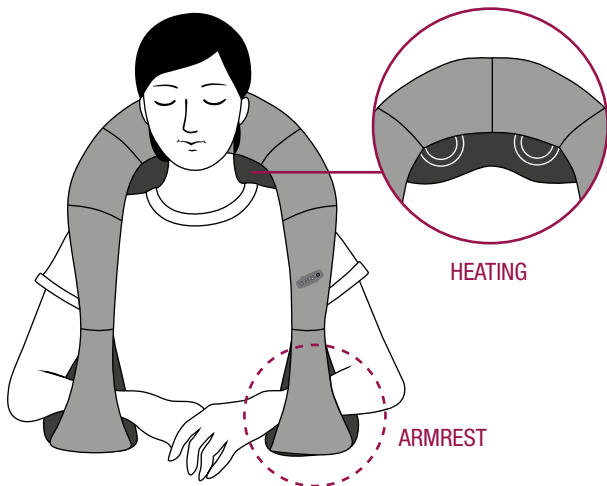 **PRODUCT FEATURES**

- The new concept of kneading massage, matching the human body shape of the neck and shoulder.
- The massage head is specially designed according to the neck and shoulder to give you the most comfortable experience.
- U shaped, arm in arm design, strong or gentle type; strength is adjustable.
- Warm function, promotes the blood circulation effectively.
- Simple switch operation.
- Polyurethane high quality material holster.
- Using the safety of the overheating protection device.

USE AND FUNCTION OF BUTTONS

Power/Heat Button

Press once, a drop of sound occurs, heat, shiatsu and tapping functions start. One press of the power button for less than 3 seconds stops the heat function; press twice and the heat function starts. Press the power button for more than 3 seconds to turn off all functions.

Shiatsu Button

When the power is on, the low level speed starts. Press the shiatsu button once (there are 2 drops of sound) then transfer to high level, the shiatsu speed will turn faster, at a speed of about 32 rounds per minute. Press the shiatsu button twice, there is a drop of sound, then transfer to low level speed. Press the shiatsu button for more than 3 seconds to turn off the shiatsu function.

Tapping Button

When the power is on, in default mode 1, there are 9 tapping modes. When transferring to another mode, a drop of sound occurs each time. When returning to the first mode, 2 drops of sound occurs. When switching into another mode, the tapping intensity will transfer to the lowest level. Press the tapping button for more than 3 seconds to turn off the tapping function.

Tapping Intensity

Press this button to switch the tapping intensity. When the power is on, the tapping intensity is at level 1. Press the tapping intensity button for more than 1 second, the intensity will transfer to the next higher level until it reaches level 9. When the intensity is at level 9, 2 drops of sound occurs. Press the button for more than 3 seconds, the intensity level will decrease each time. A drop of sound occurs each time the intensity decreases a level. When the intensity decreases to level 1, 2 drops of sound occurs.

You can use the massage device for the whole body as shown in the picture:

WAIST MASSAGE

ABDOMINAL MASSAGE

LEG MASSAGE

FOOT MASSAGE

CLEANING & MAINTENANCE

Attention: before cleaning the device, make sure that it is switched off and that the AC adapter is disconnected from the wall socket!

- After using, please use a wet cloth to clean up dust and dirt on the product.
- Please do not use corrosive cleaner, such as benzene or thinner to clean this product.
- Please do not use a wet cloth to clean the control panel, as this can cause damage. Only use a dry cloth to wipe the control panel.
- Products not used for long periods of time, should please be wipe cleaned, loaded into the case and placed in a ventilated and dry place, in order to avoid affecting the use of product due to damp.

! TROUBLESHOOTING

- If you have problems not mentioned in the list, you should try following steps.
- Turning off the power of the product, unplug the power plug from the socket
- Keep this product closed 60 minutes.
- Insert the plug into the power socket and open the switch.
- If the problem persists, please contact our customer service center.

Fault	reason	Handling
the product can't work	the plug is unplugged	plug into the socket safely
	power button is not on	turn the power on
	this product has just finished a 15 minute massage and automatically switches the machine to standby mode	if you want to massage again, please press the button
	product exposed to an overload (for example, someone standing on the product, or have a weight placed on top of the product)	unplug the product power plug, remove the excess weight. Please plug the product back in again and turn it on
Product heat	you already make the product more groups of massage 15 minutes	remove the power plug and let the product cool down. Unplug the power plug of the product after using it
The power supply cord damaged, the product smoke	cool	stop using immediately and contact the seller
The thermostat heating disconnect protection	internal explosion	

DISPOSAL

Please dispose the package accordingly. Add carton and cardboard to the waste paper and foils in the recyclables collection.

Applicable in the EU and other European countries with systems for separate collection of recyclable materials.

Do not dispose used appliances in household waste. Every consumer is legally obliged to hand in old appliances separately from household waste (e. g. at a collection point). This is the only way to ensure that old appliances are recycled properly and do not have a negative impact on the environment.

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At the following link
you will find the declaration of conformity:

www.relaxxnow.de/manuals/pdf/relaxxnow_rxm300_conformity.pdf

Distributed by:

InnovaMaxx GmbH

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service@innovamaxx.de

 **+49 30 220 663 551**

www.relaxxnow.de